
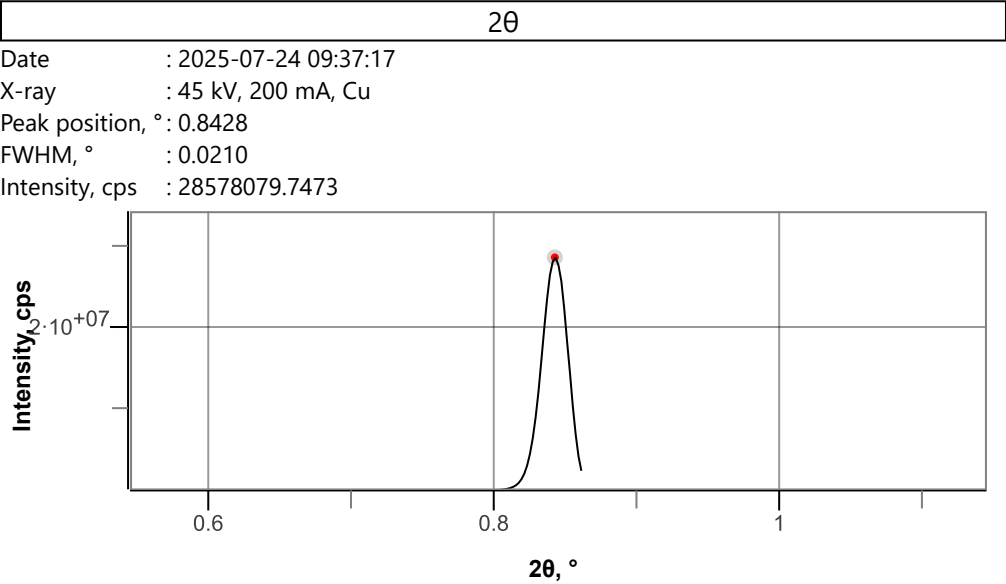
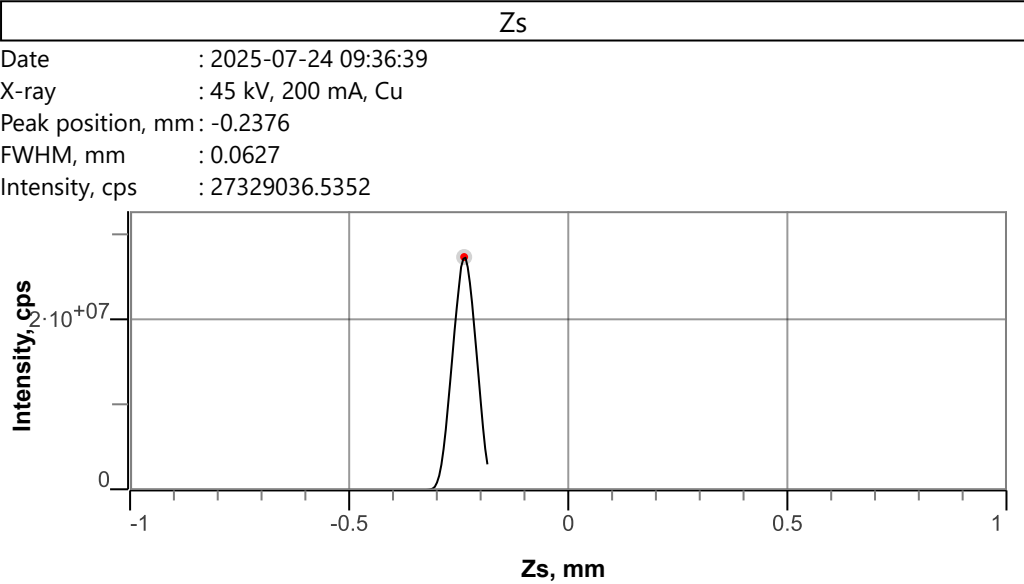
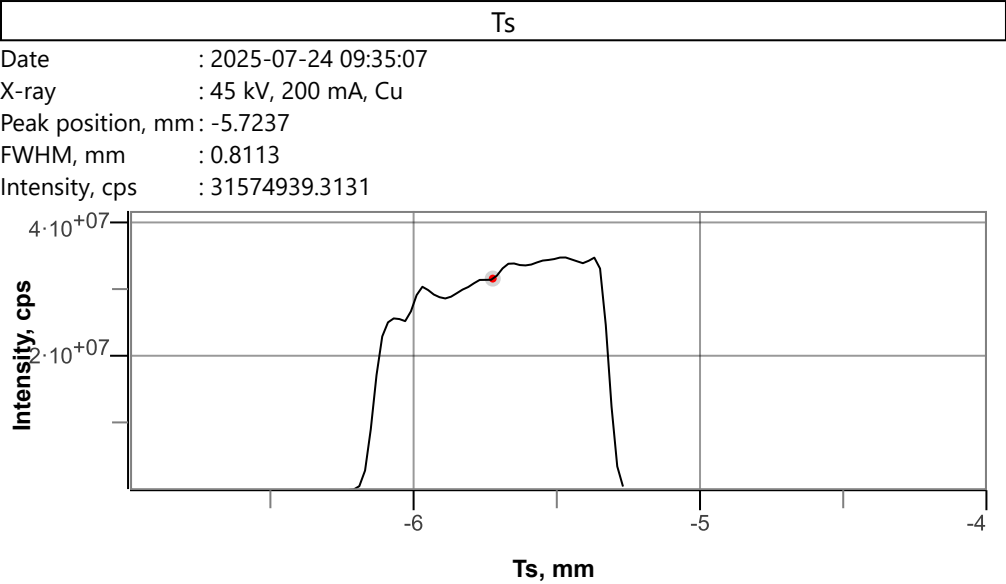


Optics Alignment Results

X-ray generator : 9kW Rotating anode(45kV 200mA ID:1703)
Target product id : Cu Fine focus (A-26-Cu)
Goniometer : Inplane Goniometer with encoder
Attachment base : $\chi\phi Z$ attachment
Attachment head : Standard attachment head
Incident primary : Standard
CBO : CBO
Selection slit : PB
Incident optics unit : Ge(220)x2
Incident Soller slit : Soller slit open (short)
Incident slit : Integrated incident slit box
Length-limiting slit : 10 mm
Receiving slit #1 : Receiving slit box #2
Filter 1 : None
Receiving optics unit #1 : PSA open
Receiving Soller slit : Soller slit 2.5°
Receiving slit #2 : Virtual receiving slit box
Filter 2 : None
Receiving attenuator : No unit
Detector : HyPix-3000 (horizontal)
Detector monochromator: None
Detector slit : None


Optics attribute : PB-Ge(220)x2
Optics alignment name: User defined settings

Axis	Position	Offset
2 θ		0.8428
ω		0.7805
Ts	-5.7236875	
Zs	-0.2375000	
Zr	0.1471875	
2 $\theta\chi$		0.0000
M	0	
ωm	0	

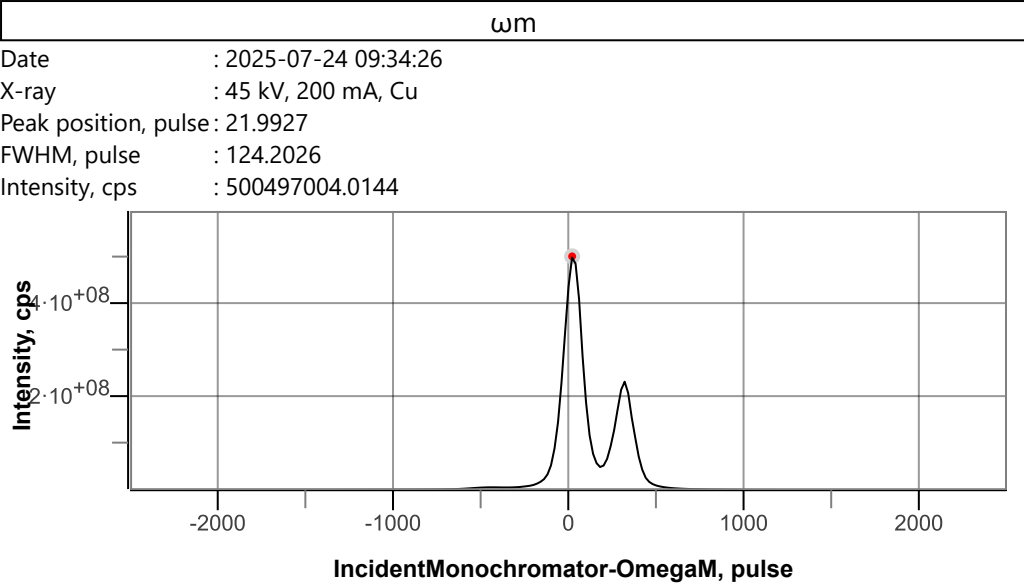

Zr

Date : 2025-07-24 09:37:56
X-ray : 45 kV, 200 mA, Cu
Peak position, mm : 0.1471
FWHM, mm : 0.0653
Intensity, cps : 23491707.2706

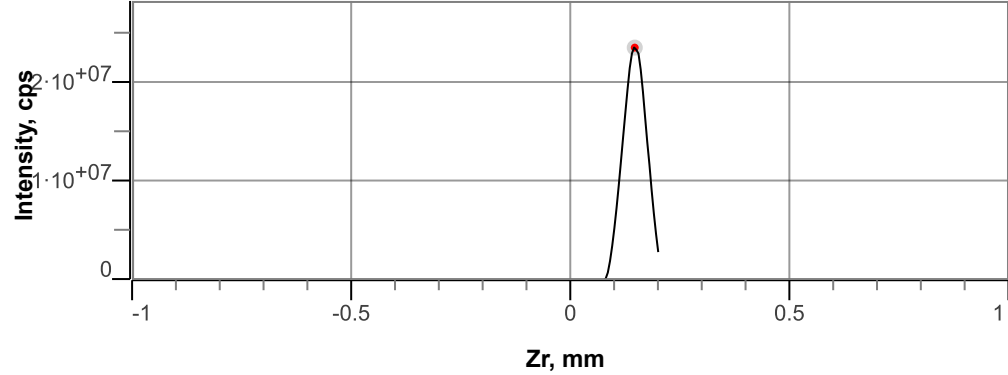

ω

Date : 2025-07-24 09:38:38
X-ray : 45 kV, 200 mA, Cu
Peak position, ° : 0.7805
FWHM, ° : 0.1723
Intensity, cps : 21588403.7690

